

Mini Fascia Gun

C503-8

Operating Instructions

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the manual, it is a normal update of the product, and our company will not inform you otherwise. The specific product function and color shall prevail in kind.

CONTENTS

Important Safety Warning	02
Name and function of components	04
Method of usage	05
Troubleshooting	07
Operate Instruction	08
Transport and storage	08
Name and content of harmful substances	09
Product specification	10

Important Safety Warning



1. Safety Warning

- This instrument is not suitable for people (including children) who are physically disabled, sensorily or neurologically defective, or lacking experience and common sense in using such instruments, unless they are supervised and guided by those responsible for their safety on how to use the instrument. Children must be monitored to ensure that they will not play with this instrument.
- Please follow this manual to use this instrument.
- Please use a charging cable that matches this instrument or use a charging cable with a charging current of 2A or more.
- It is forbidden to splash liquid on this product to avoid accidents caused by short circuits.
- Please pull out the power plug after use or before cleaning to avoid damage to the equipment or injury personnel.
- Do not insert or pull out the USB charging cable with wet hands.
- Do not apply an excessive force to the USB charging cable when inserting or pulling it out.
- Do not damage the wires or tamper with the circuit of this instrument.
- The massage time for the same part should not exceed 60 seconds. For those with sensitive skin, cover the massage part with a towel or put on clean and dry clothes before use.
- Do not cover the air vent of this product.
- It is forbidden to throw the product into the fire or put it near a fire source.
- Any modification of this product without authorization is prohibited.
- Keep your fingers, hair or other parts of the body away from the massage shaft and heat dissipation holes so as to avoid being caught.
- Do not apply great pressure on, knock or smash the machine during use.
- Do not use accessories and fittings that are not recommended.
- If you feel any abnormality of this instrument during use, stop using it immediately and consult your local supplier.

Important Safety Warning





2. Working Environment

- Do not use this product in a highly humid environment such as a bathroom.
- Do not use this product immediately after there has been a sharp change in the environmental temperature.
- Do not use this product in an environment where there are lots of dust or corrosive gases.

3. People That Should Not Use This Product

- Women that are pregnant or menstruating.
- People with contagious skin conditions.
- Patients with osteoporosis.
- Patients with heart diseases and with electronic medical instruments like a heart pacemaker.
- People in a high fever.
- Those who are injured or suffering from anybody surface disease.
- Children under 14-year old or people mentally unbalanced without being supervised are not allowed to use this machine.
- Those who need to rest according to a doctor's advice or those who do not feel well.
- Those who are wet.
- Those who are under the influence of alcohol or feeling unwell.
- Please according to your doctor's advice if you have had an operation recently.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.


Name and function of components

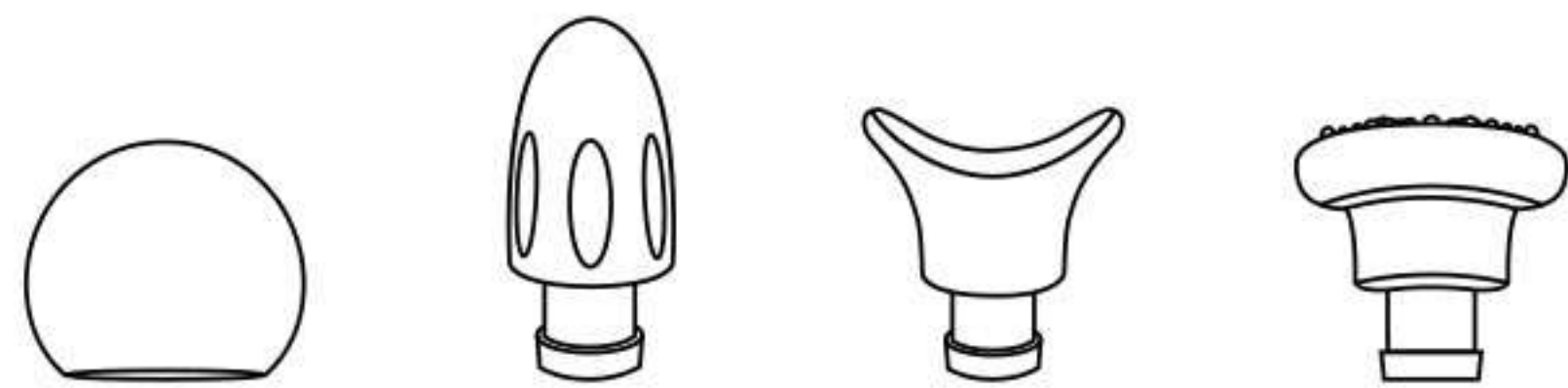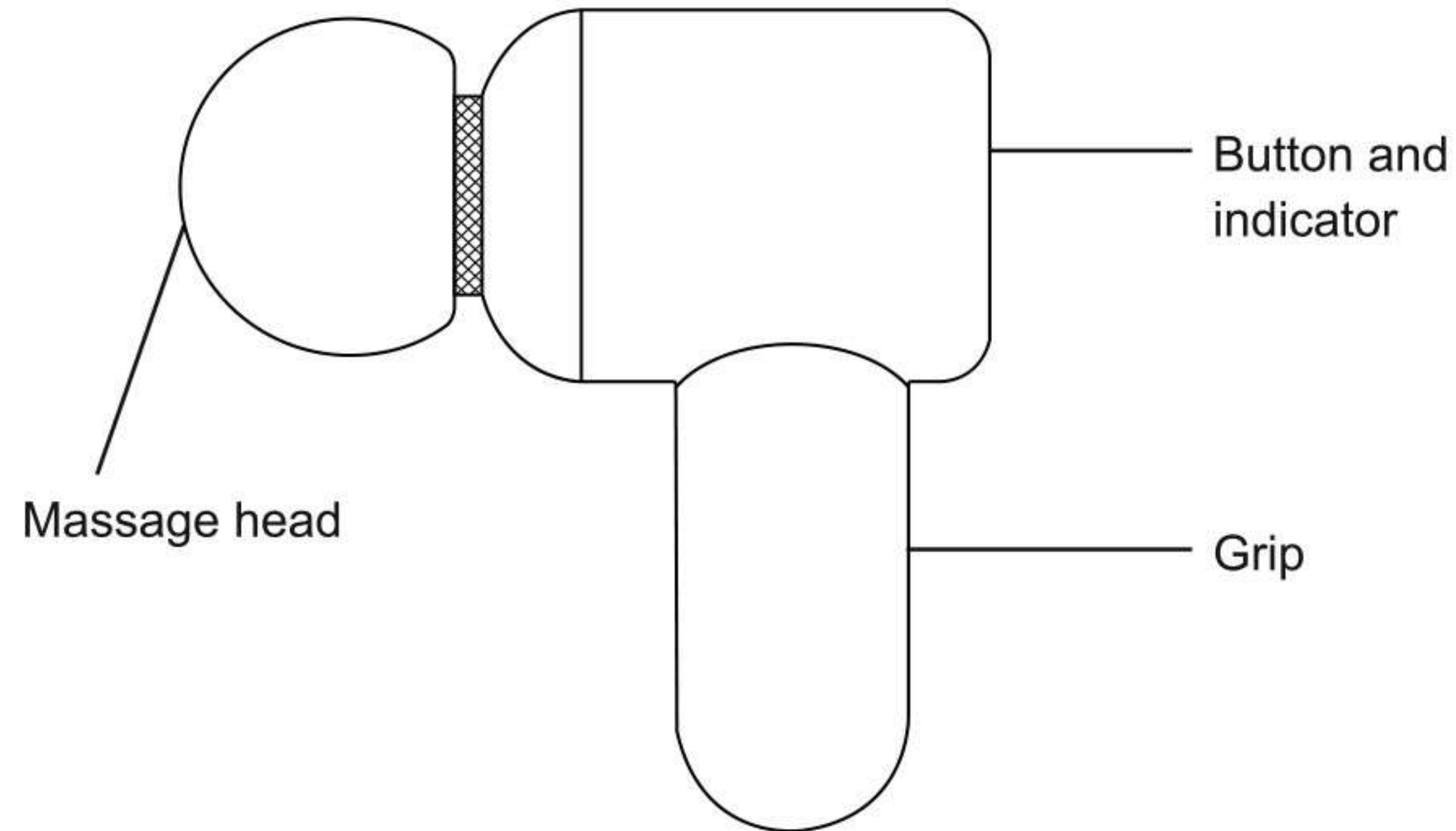

- Spherical head: suitable for massaging big muscle groups such as arms, back, buttocks, thighs and calves.
- Cylindrical head: suitable for impacting deep tissues such as meridians, palms and soles.
- U-shaped head: suitable for massaging the neck, spine and Achilles tendon.
- Flat head: suitable for relaxation and shaping of various muscle parts.

SAVE THESE INSTRUCTIONS

Method of usage

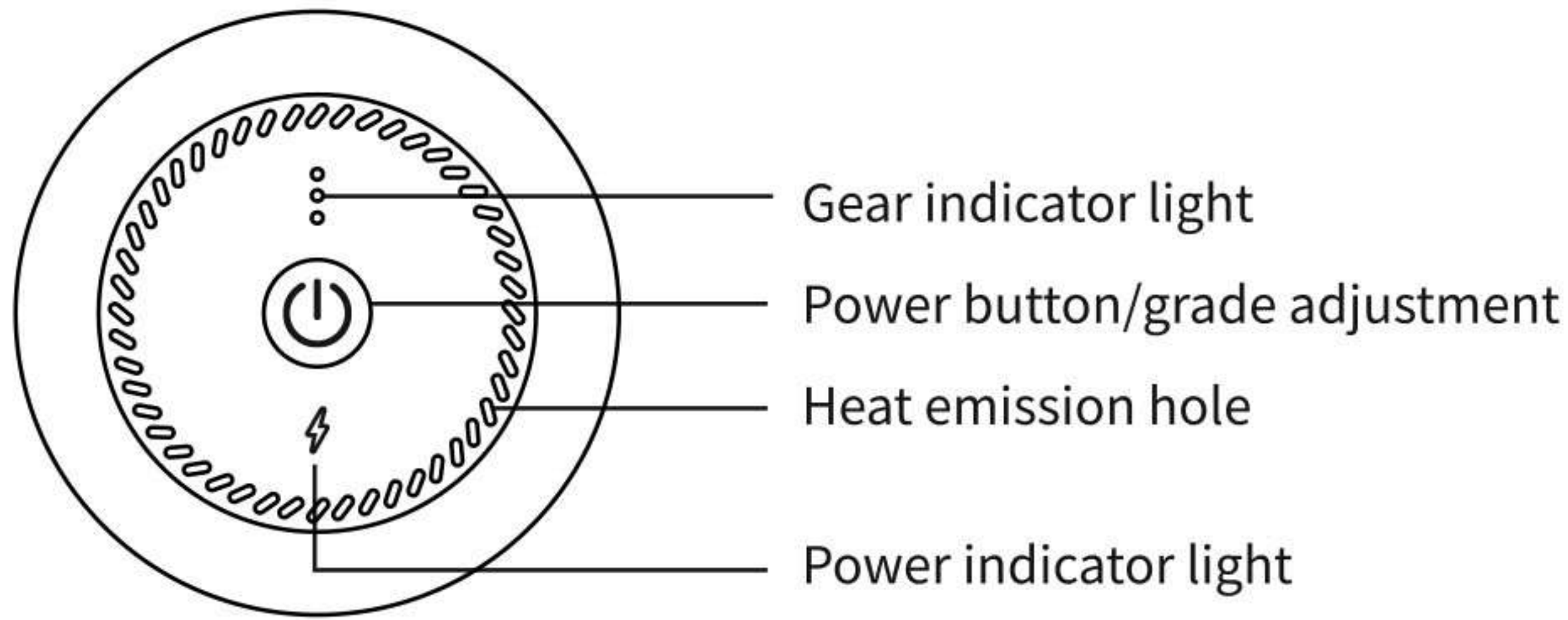

Instructions for use:

- The raised rubber button in the middle of the control panel is the power button (speed control button). Press and hold it to turn the instrument on or off.

Turn on the massager, it starts with level 1. Press the power button continuously to change the speed from level 1 to level 6 and to level 1 again.

Press the button quickly for two times, and all the three lights flashing. The massager is in intelligent mode, it adjusts the speed automatically according to the strength you press on it. When the massager is turned off, short press to show the remaining battery capacity.

The indicator flashes slowly in red, indicating that the instrument is being charged, while the indicator is on in green, indicating that the instrument is fully charged.

•

- 0-5% The indicator flashes quickly in red (the instrument is prohibited to be turned on and should be charged as soon as possible);
- 5-30% battery level, the indicator is red; 30-60% battery level, the indicator is yellow; and 60-100% battery level, the indicator is green.

Method of usage

- In order to avoid massaging the same part for too long time and damaging such muscles and bones, the machine will enable the intelligent protection and stop working after 10 minutes of continuous work. To you use the machine again, you need to press the power button.
- The charging port is at the bottom of the button indicator panel.

Through strong and beneficial high-frequency vibration stimulation, this product can penetrate deep into the texture, penetrate the deep fascia, quickly disintegrate the lactic acid accumulated after exercise and relax tight and rigid muscle tissues. While combing knotted muscles and preventing muscle adhesion, it completely remove soreness and relieves fatigue and discomfort.



Muscle groups of shoulder
15s~20s x 3



Muscle groups of chest
30s



Muscle groups of back
30s~1min x 3



Muscle groups of buttocks
30s~1min x 4



Muscle groups of legs
20s~30s x 3



Muscle groups of feet
20s~30s x 3

- This product is only used as a physical therapy and health care product, not as a medical appliance.

Method of usage

Contraindications:

- Without the permission of a doctor, this product must not be used in any of the following situations:
- If you are pregnant, suffer from complications of diabetes (such as neuropathy or retinal damage), epilepsy or migraine, herniated disc, spondylolisthesis, spondylolysis or spondylosis, wear a pacemaker, recently had joint replacement, IUD, metal needle or steel plate or other surgery and have any other problems related to your health.

Troubleshooting

Basic Fault	Reason	Fault solution
Can not turn on	Power off	Charge the product if is power off
Can not charge	Adaptor is good connect with wire or not	Right connect the wire and adaptor
	Wire or USB plug is damaged or not?	If damaged, contact the manufacturer
	Is the adapter used in compliance with usage requirements?	Use the charger specified by the manufacturer
The product can be started, but the function can not be started	Check the remaining charge.Is the charge weak?	recharged

If a fault occurs when operating this product, please check according to the following methods, if cannot find the problem or the product still cannot operate, please shut off the power and contact the service center to repair.

Operate Instruction

Product battery dismantling steps:

- 1.Cut off the external power supply and remove the product shell;
- 2.Cut the two power cords connecting the battery in turn;
- 3.Remove the battery, insulation glue stick to the power cord is broken respectively;
- 4.Please place the discarded batteries at a fixed recycling point;
- 5.Please follow local environmental laws and guidelines to discard and care for the environment.

Note:

- The battery must be removed from the appliance before discarding it;
- The appliance must be powered off when the battery is removed;
- Batteries should be disposed of safely;
- Do not dismantle it for non-professionals.

Transport and storage

1.Temperature:

Keep the whole machine in a dry and cool place.

2.The relative humidity: no more than 80%, no corrosive gas and Well-ventilated room.

3.Transport: This product can bare simple shock should be avoid to raining, snowing, raining splash, inversion and collision. If the product storage more than six months, it should be check and c harge half hour to maintenance.

The room should be well ventilated , avoid intense sun exposure and corrosive gas.

Name and content of harmful substances

Declaration of Toxic and Harmful Substances in Products

Component name	Toxic and harmful substances or elements					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated diphenyl ethers (PBDE)
Main body of massage device	○	○	○	○	○	○
Data wire	○	○	○	○	○	○

The table is created according to SJ/T11364-2014.

- :Indicates that the content of the harmful substance in all homogeneous materials of the component is below the limit specified in GB/T 26572-2011.
- : Indicates that the content of the harmful substance in at least one of the homogeneous materials of the component exceeds the limit specified in GB/T 26572-2011.



This product is free of harmful substances, is green and environmentally friendly. It can be recycled and should not be discarded at will.

Product specification

Model: C503-8

Name: Mini Fascia Gun

Rated time: 10min

Rated power: 20W

Rated voltage: 5V $\overline{=}$ 2A

Charging interface: Type-C

Lithium battery capacity: 7.4V 2200mAh

Charging time: \approx 3h

Battery life: \approx 4h

Product safeguard construction: 